



# REVOLUTIONIZING PATIENT REPOSITIONING!!



**PRESSURE RELIEVING &  
REPOSITIONING SYSTEM**



# PATIENT COMFORT AND SAFETY

**"Say goodbye to manual patient repositioning!"**

The automated system repositions patients every 30 minutes effortlessly!!"

Say goodbye to the need for repositioning every two hours; EVA TILT is here to take on that responsibility, freeing up valuable time for nurses and transforming the patient experience from pain to joy.

Our innovative tilt system places patient comfort as our top priority, along with recognizing the significant time invested by caregivers in positioning patients. With EVA TILT, ensuring optimal patient positioning is just a touch of a button away.

EVA TILT goes beyond the standard care by automatically adjusting pressure points every five minutes and repositioning the patient every 30 minutes - exceeding the international recommendation of repositioning every two hours.

SAVING  
NURSE'S  
TIME



EVA TILT It's a Game-changer in Patient Care.

What sets EVA TILT apart is its variable low-pressure mode, providing an advanced support surface that surpasses traditional alternating pressure mattresses. This feature ensures a higher level of pressure relief and redistribution, elevating the standard of care for patients.

Bid farewell to the reliance on traditional alternating pressure mattress systems; embrace the future with EVA TILT's advanced repositioning support surface.

Dynamic mattress replacement suitable for users at "High Risk" of Pressure Injury.



## EASY PATIENT RE-POSITIONING



BOTH ALTERNATING AND TILTING WORKS IN TANDEM



## A SMART SOLUTION FOR PATIENT CARE.

Our EVA TILT mattress automatically shifts pressure points every five minutes. With a gentle tilt of up to 30 degrees in every 30 minutes, which ensure optimal pressure relief and compliant patient repositioning.



Adjusting pressure points every 5 minutes




Repositioning every 30 minutes



Tilts upto 30° angle



Elevate patient care, save time and transition from pain to joy with EVA TILT!

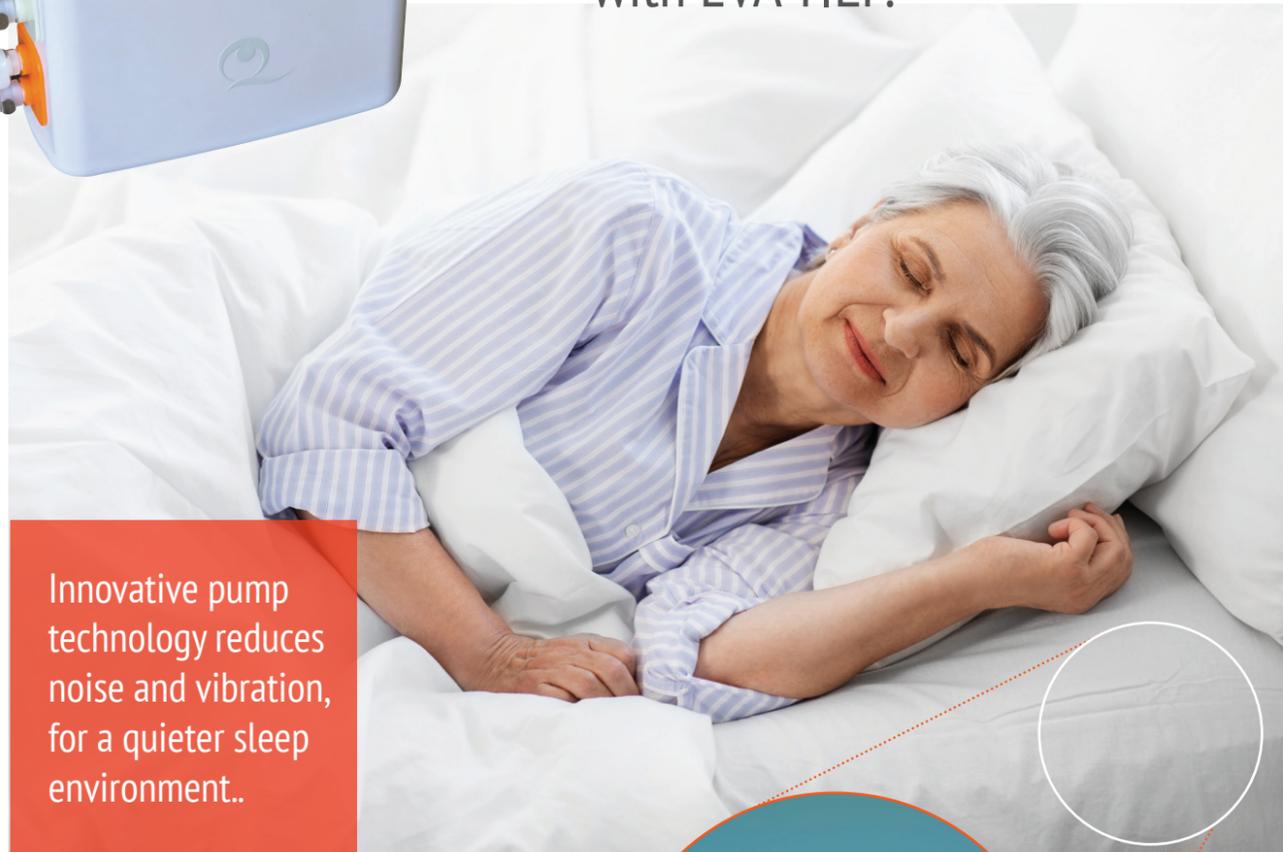

Innovative pump technology reduces noise and vibration, for a quieter sleep environment..



**LOW FRICTION PU COVER**

Low friction cover which is made of polyurethane cover reduces shear force.



**WATER PROOF COVER**

The cover is washable making it most hygienic in managing patient affairs. Fully water repellent and easy to clean is one of the key features in preventing infections.

**INDIVIDUALLY REPLACEABLE CELLS**

The EVA® offers a replaceable cell which can be quickly changed to keep the mattress's life extended to many years. The cell construction makes it possible to completely remove the individual cells for cleaning and decontamination.



**USER INTERFACE**

EVA® ensures simplicity and ease of use. The dial knob and color coding offers easy set up and patient comfort.



**CPR FUNCTION**

Fast Access to Immediate CPR A simple pull of the CPR tag provides firm support for commencing chest compressions..



**NON-INTERCHANGEABLE QUICK CONNECTOR**

A unique non-interchangeable quick connector with CPR function makes the connection between mattress and control unit simple and easy to use.



**DETACHABLE POWER CODE OPTION**



**SOFT FRICTION HANGERS**

The soft pump hooks fasten effectively to a wide range of beds.



**WELL PROTECTED ZIPPER**

The Cover is attached with concealed two-way zips to reduce the ingress of fluids.

**ADVANCED DIGITAL CONTROL UNIT**





### SPECIFICATION FOR EVA TILT

Model	CXIP40
Input voltage	120-240V/50-60Hz
Mode of operation	VLP, CLP, Maxfirm
VLP cycle time	12-15min
Tilt cycle time	90Min
No. of repositioning in tilt cycle	4
Power consumption	0,010Kw/h
Weight	3kg
Working temperature	10-50
Humidity	Less than 95%
Electrical classification	Class 3 type BF
Applied parts	CXIM mattress

#### MATTRESS

Model	CXIM40
Cell type	Nylon TPU
Outer cover	Polyurethane coated fabric
Dimension	90*200*10
Maximum user weight	150kg



### SPECIFICATION FOR EVA TILT PLUS

Model	CXIP42
Input voltage	120-240V/50-60Hz
Mode of operation	Auto, Pulsating, CLP, Maxfirm
Cycle time	12-15min
Tilt cycle time	90Min
No. of repositioning in tilt cycle	4
Power consumption	0,010Kw/h
Weight	3kg
Working temperature	10-50
Humidity	Less than 95%
Electrical classification	Class 3 type BF
Applied parts	CXIM mattress

#### MATTRESS

Model	CXIM42
Cell type	Nylon TPU
Outer cover	Polyurethane coated fabric
Dimension	90*200*20
Maximum user weight	220kg

AVAILABLE IN BOTH OVERLAY AND  
MATTRESS REPLACEMENT MODELS

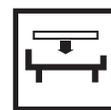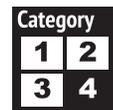

Always read the user manual prior to use.



CEUMED MEDICAL DEVICES LLP

SDF-2, KSIDC Investment Zone, Kanalpirivu, Palakkad,  
Kerala-678621, INDIA @ info@ceumed.com

www.ceumed.com

